


## Leksell® Multi Purpose Stereotactic Arc Kit - with Counter Scale

Article No : 1008174

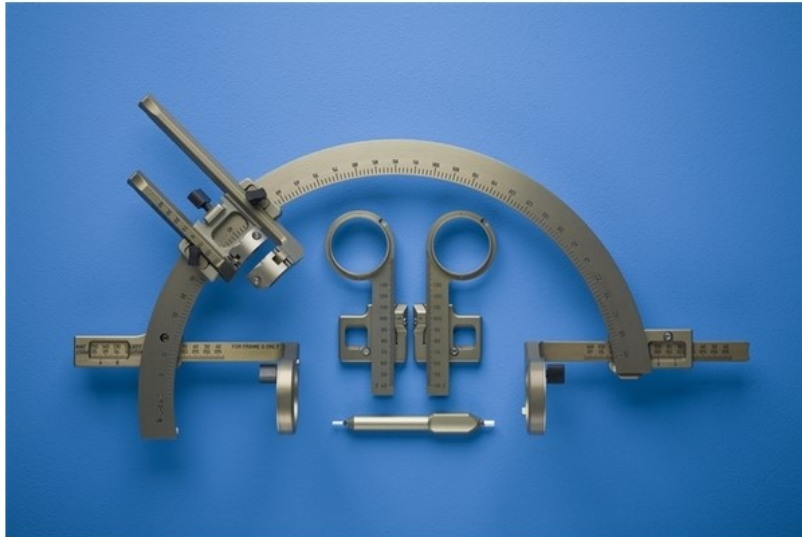

### All-around surgical approach with ensured accuracy.

The x-coordinate is set on both axes before mounting, which assures fast and correct positioning when fitting the arc to the frame.

Leksell® Multi Purpose Stereotactic Arc with counter scale is designed to allow the surgeon to select any desired trajectory to the predetermined target, regardless of the patient's position, when using Leksell Coordinate Frame G.

Thank's to the counterscale the accuracy is further ensured when mounting the arc.

Leksell® Multi Purpose Stereotactic Arc with counter scale can easily be mounted onto the frame in either A-P or lateral orientation. Lateral orientation is most commonly used, but the A-P orientation is especially useful for the implantation of depth electrodes in lateral targets.

The unique option to attach several sliding carriers on the arc in combination with the center-of-arc principle allows the surgeon to simultaneously reach and/or observe any selected target through two or more trajectories.

### Consisting of

- Leksell® Multi Purpose Stereotactic Arc with counter scale
- Slide carrier to Leksell® Multi Purpose Stereotactic Arc
- Arc Support with Slide, Right
- Arc Support with Slide, Left
- Instrument Screw Driver, Double
- Knob Spare Part Kit
- MPS Arc Spare Part Kit
- Sterilizing Tray for Multi Purpose Arc
- Quick Reference Guide, Leksell® Multi Purpose Stereotactic Arc with Counter Scale
- Stainless Steel Screw Kit M4x8
- Nylon Screw Kit M3x4