


MR-Adapter for GE Split Headcoil

Article No : 1020954

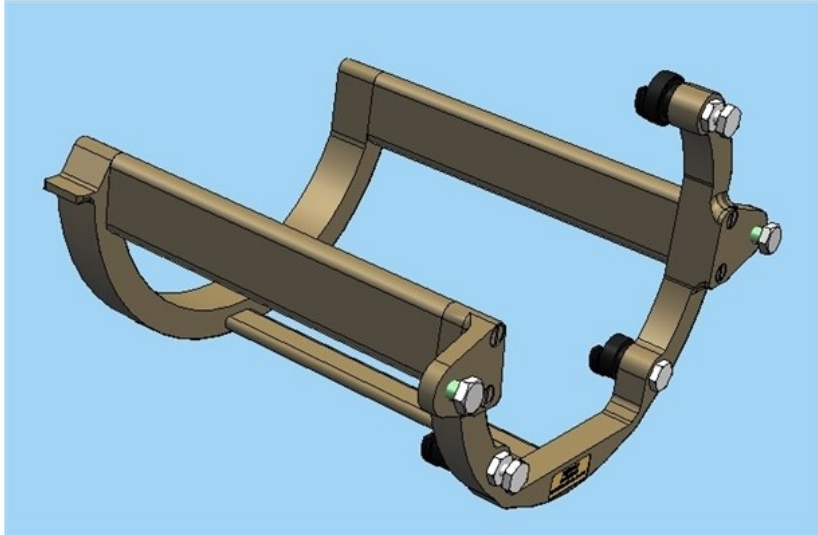

The adapter is designed for GE 450/750 Split Head coil and GE HD Split head coil with inner diameter of 281 mm and a length of 424 and 414 mm respectively.

Each adapter comprises a base which attaches to the MR scanner table and provides a rigid platform for the frame support. The MR adapter is constructed from a nonmagnetic material. The exact appearance of the adapter and the means by which it is attached to the scanner table vary according to the type of scanner. A range of MR adapters to suit different designs of MR scanners is available from Elekta.

This part replaces MR Adapter GE HD Split Head Coil with art no 917780.

For further details on which adapter you require for your scanner, please contact Elekta.