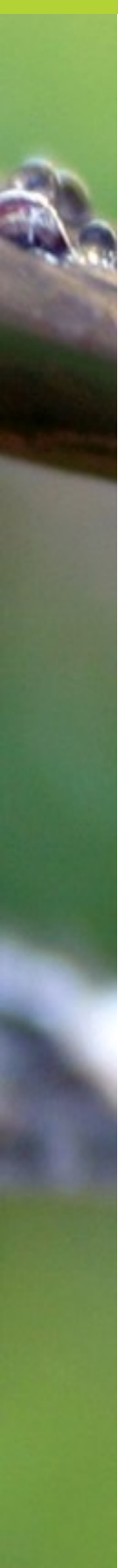


Table Fixation for X-Ray, Flat 228 mm

Article No : 407010

A vertical strip on the left side of the page shows a close-up of a patient's head with a coordinate frame attached. The frame consists of several small, dark, circular markers arranged in a grid pattern on the patient's scalp.

Designed for flat table head end, width 228 mm. For tables with a thickness of minimum 30 mm and 45 degree angle on the table sides. X-Ray Table Fixation, Table width 228 mm The X-Ray Table Fixation is fitted to the table of the X-Ray machine. The table fixation is designed to hold the coordinate frame, indicator and the skull of the patient firmly and to ensure that position and orientation can be reproduced in subsequent scans.

A range of X-Ray table fixations to suit different designs of X-Ray tables are available from Elekta.

For further details on which table fixation you require for your table, please contact Elekta.