


## Sedan Side-Cutting Biopsy Kit 2 - 5 mm

Article No : 1522155

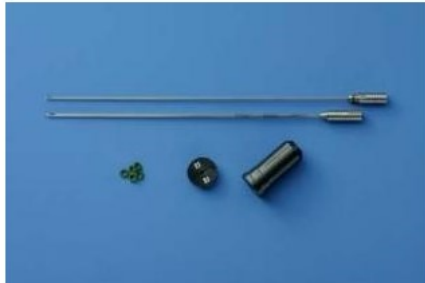

New biopsy needle with a 5 mm needle window. The vacuum and side-cutting principle of Sedan has been adapted to Leksell Stereotactic System® and Leksell® Vantage™ Stereotactic System in this biopsy needle. Sterilizing tray and Stop and Guide Inserts with a diameter of 2.1 mm are included in the kit.

### Accurate sampling for intact specimen

Sedan Biopsy Needle Kit 2 is designed for use with Leksell Stereotactic System® or Leksell® Vantage™ Stereotactic System. The working length is the Leksell Stereotactic System® standard (190 mm).

Sedan Biopsy Needle Kit 2 comprise a set of neurosurgical instruments designed to allow multiple biopsies from one or more targets in one trajectory. The sample is obtained by applying vacuum. The window in Kit 2 is only 5 mm in length, and so allows a smaller specimen to be obtained. The 5 mm window size permits that a series of specimen to be taken.

The kit also contains 5 O-Rings 3.1 x 1.6 mm Viton.

This part replaces A2430-02.

Manufacturer: Elekta Instrument AB

### Consisting of

- Sedan Side-Cutting Biopsy Needle - 5 mm
- Guide with 2,1 mm hole
- Stop with 2,1 mm hole
- Sterilizing Tray / Small