


## Leksell GammaPlan® -Inverse Planning

Article No : 1016650

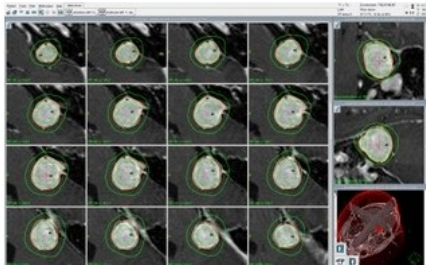

This software add-on enables fast automatic planning. It consists of two key independent functions –filling and optimization – which affect the critical activities of dose distribution and plan optimization. With “filling”, the target is automatically filled with shots. The optimization function is used to improve a preliminary plan. The optimization process is controlled by prioritizing four parameters: target coverage, selectivity, dose gradient and treatment time. It is also possible to change the number and location of shots at any time during the procedure. By making the optimization very fast, an efficient iterative workflow is supported, where the plan is developed through adjustment of the inverse planning settings, continued optimization and possibly manual planning.